

SYNAPSE

DENTAL PAIN ERASER

HOW TO USE AT HOME



TURN ON

Press button to turn on light.

BRACES PAIN?

Move tips on gums (like painting a pole) for ~20-30 sec per tooth.



TOOTH PAIN?

Move tips ~2 min on gums (like painting a pole - see above) and 10 sec on tooth.

MOUTH SORES?

Full circle AROUND sore for ~30 sec.

NO tips ON sore.



TURN OFF

Press button 2 sec.
Light flashes then off.
Auto shut off in 5 min.

CLEAN

Rinse and dry tips.
Replace cap.



LEARN
MORE



TIPS FOR USING



USE UNTIL PAIN GOES AWAY
20 sec up to 5 min.

BOTH TIPS MUST TOUCH GUMS OR TOOTH!

Use with a mirror or while being still to make sure (watching TV, reading, etc.)



DON'T SHARE
Use only in your own mouth.

SALIVA HELPS

Too dry or too much saliva? Close lips and swallow.



PREVENT PAIN
Use BEFORE your dental visit, too!

Once you've used the Dental Pain Eraser, it will change your dental experience forever!



SYNAPSE DENTAL PAIN ERASER



REWRITING THE DENTAL EXPERIENCE

Prevent and relieve oral pain

Take one home today!


PATIENT GUIDE TO HOW IT WORKS

AAO Ortho Innovator Award
& Best New Hygiene Product



DELIVERY OF IMMEDIATE DENTAL PAIN RELIEF →

Teeth, cheeks, lips, tongue and gums contain nerve endings that are highly sensitive for pain, temperature, air and pressure. In fact, a large part of our brain is used for sensory perception in our mouth.

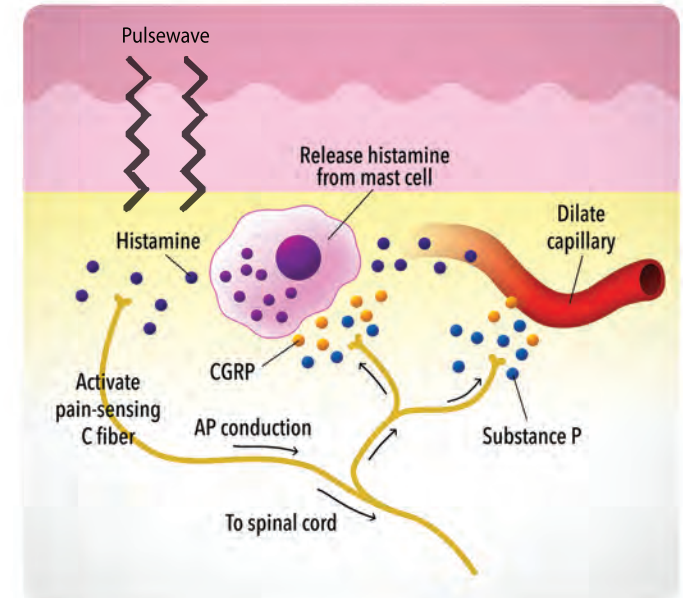

THE IDEAL PAIN MANAGEMENT TOOL →



Erase Dental Pain

Million of patients avoid the dentist because they are afraid. With the Dental Pain Eraser, you never have to be anxious about dental care again!

Apply on root surface and on tooth surface for up to 45 seconds per tooth on average.



Erase Canker Sore Pain

1 in 5 people experience canker sores routinely. Immediate abrasion relief from canker sores, aligners, braces, or dentures!

Apply around the canker sore or abrasion for 20 seconds, not touching the inflamed area.

DELIVERY OF IMMEDIATE DENTAL PAIN RELIEF

(Continued)

Advanced Pulsewave™ technology sends gentle electric pulses to calm the nerves in the mouth and provide long-lasting pain relief. Can also be used for **pain prevention**.

THE IDEAL PAIN MANAGEMENT TOOL

(Continued)



Erase Ortho Pain

Proven immediate pain relief during orthodontic treatment. Pain will never be an obstacle to a beautiful straight smile again!

Apply up and down the root surface for 10 seconds per tooth with an average of 2 minutes for the whole mouth.

IMPORTANT SAFETY INFORMATION

The Synapse Dental Pain Eraser is not suitable for everyone. Please do not use this device if you have a pacemaker, implanted defibrillator, or any other implanted electronic or metal device, as it may interfere with their function.

Avoid use if you have a history of cerebrovascular conditions, heart conditions, epilepsy, or are prone to seizures. This device should also not be used for unexplained or undiagnosed pain.

If you have any of these conditions or are unsure whether this device is appropriate for you, please consult your doctor or orthodontic provider before use.

"Highly Recommended by Clinicians for Their Patients" by Clinicians Report®, the Top Independent Dental Tech Evaluator



“My daughter is vocal when it comes to pain. After using the Dental Pain Eraser, she wiggled her loose tooth and it didn't hurt anymore. The Dental Pain Eraser did the trick! I love using it as a mom and orthodontist.

Laura Bonner, DDS, MS
Diplomate, American Board of Orthodontics



Synapse Dental, LLC

Synapse Dental is a dental solution provider committed to making pain-free dentistry a reality for dental professionals and patient.



dentalpaineraser.com

Rev B, 05-23